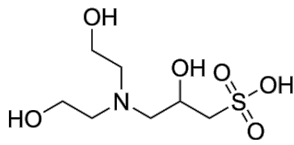


DIPSO

货号: **B26940**

产品信息

生物活性	DIPSO is a biological zwitterionic buffer with the useful pH range from 7.0 to 8.2. DIPSO can interfere meiotic regulation in mouse oocytes. DIPSO also has surfactant activity at 10 mM.
CAS	68399-80-4
中文名称	
分子量	243.28
体外研究	
体内研究	
形式	Solid
运输条件	Room temperature in continental US; may vary elsewhere.
保存条件	

溶解性	In Vitro: H₂O : 166.67 mg/mL (685.10 mM; Need ultrasonic) 配制储备液 <table><tr><th>浓度</th><th>溶剂</th><th>体积</th><th>质量</th></tr><tr><td></td><td></td><td>1 mg</td><td>5 mg</td><td>10 mg</td></tr><tr><td>1 mM</td><td>4.1105 mL</td><td>20.5525 mL</td><td>41.1049 mL</td><td></td></tr><tr><td>5 mM</td><td>0.8221 mL</td><td>4.1105 mL</td><td>8.2210 mL</td><td></td></tr><tr><td>10 mM</td><td>0.4110 mL</td><td>2.0552 mL</td><td>4.1105 mL</td><td></td></tr></table> * 请根据产品在不同溶剂中的溶解度选择合适的溶剂配制储备液；一旦配成溶液，请分装保存，避免反复冻融造成的产品失效。 储备液的保存方式和期限：-80℃, 6 months; -20℃, 1 month。-80℃ 储存时，请在 6 个月内使用，-20℃ 储存时，请在 1 个月内使用。 In Vivo: 请根据您的实验动物和给药方式选择适当的溶解方案。以下溶解方案都请先按照 In Vitro 方式配制澄清的储备液，再依次添加助溶剂： ——为保证实验结果的可靠性，澄清的储备液可以根据储存条件，适当保存；体内实验的工作液，建议您现用现配，当天使用； 以下溶剂前显示的百分比是指该溶剂在您配制终溶液中的体积占比；如在配制过程中出现沉淀、析出现象，可以通过加热和/或超声的方式助溶 1. 请依序添加每种溶剂： PBS Solubility: 50 mg/mL (205.52 mM); Clear solution; Need ultrasonic *以上所有助溶剂都可在 MCE 网站选购。	浓度	溶剂	体积	质量			1 mg	5 mg	10 mg	1 mM	4.1105 mL	20.5525 mL	41.1049 mL		5 mM	0.8221 mL	4.1105 mL	8.2210 mL		10 mM	0.4110 mL	2.0552 mL	4.1105 mL	
浓度	溶剂	体积	质量																						
		1 mg	5 mg	10 mg																					
1 mM	4.1105 mL	20.5525 mL	41.1049 mL																						
5 mM	0.8221 mL	4.1105 mL	8.2210 mL																						
10 mM	0.4110 mL	2.0552 mL	4.1105 mL																						
纯度	≥98.0%																								